

# Arm Type Blood Pressure Monitor

Model: # BPM67B

## INSTRUCTION MANUAL

PLEASE READ THIS INSTRUCTION MANUAL COMPLETELY  
BEFORE OPERATING THIS UNIT

EN



**IP 22** RoHS REACH



MDSS GmbH  
Schiffgraben 41  
30175 Hannover Germany



AViTA Corporation  
9F., No. 78, Sec. 1, Kwang Fu Rd.,  
San Chung Dist., New Taipei City 24158, Taiwan,  
China



China, AViTA(WUJIANG)



Importer



Distributor

## Contents

Intended Use.....	1
Important Information before Use .....	2
Product Identification .....	6
Description of LCD Display .....	7
Battery Installation .....	8
Setting the Date and Time.....	9
Placement of the Pressure Sleeve.....	10
Measurement of Pulse Rate and Blood Pressure .....	12
Hypertension Indicator .....	13
Irregular Heartbeat Detector .....	14
Atrial Fibrillation (AFIB).....	14
Memory Function .....	16
Data Transferring .....	17
Error Codes .....	19
Troubleshooting .....	20
Cleaning and Disinfecting.....	21
Technical Specification .....	22
EMC Tables.....	23
Explanation of Symbols.....	25

## **Intended Use**

The product automatically measures human being's Systolic, Diastolic blood pressure and pulse rate by oscillometric method. The measurement results are displayed on the LCD. Measurement position is at human being's upper arm. The intended use of this over-the-counter device is for adults with upper arm circumference ranging from 220 mm to 420 mm (Approx.8.7 ~ 16.5 inches) and for home use. When the device detects the appearance of irregular heartbeats such as atrial or ventricular premature beats during measurement, an indicated symbol will appear with measuring readings.

This device is designed only for adults.

## **Type of Use/ Reuse**

Multiple patient multiple use

## **Intended User**

The Arm type blood pressure monitor is intended for both professional and consumer, and the patient is the intended operator.

Patient selection criteria: Handicapped persons and children are the exception, as of handicapped persons and children need assistance by another person to use the device.

## **Contra-indications**

Do not use in these cases: common arrhythmias such as atrial or ventricular premature beats, atrial fibrillation, arterial sclerosis, poor perfusion, diabetes, age, renal diseases, patient motion, trembling, shivering.

## Important Information before Use

1. Blood pressure measurements should only be interpreted by a physician or a trained health care professional who is familiar with your medical history. Through regular use of this device and recording of your measurements, you can keep your physician informed of the changes in your blood pressure.
2. Perform your measurement in a quiet place. You should be seated in a relaxed position.
3. Avoid smoking, eating, taking medication, alcohol consumption or physical activity 30 minutes prior to making a measurement. If you are exhibiting signs of stress, avoid taking your measurement until the feeling subsides.
4. Rest 15 minutes prior to taking a reading.
5. Remove any constrictive clothing or jewelry that may interfere with the cuff placement.
6. Keep the monitor stable during measurement to achieve an accurate reading. Remain still; do not talk during the measurement.
7. Record your daily blood pressure and pulse readings on a chart.
8. Take your readings at the same time each day or as recommended by your physician to get an accurate indication of change in your true blood pressure.
9. Wait a minimum of 15 minutes between readings to allow blood vessels to return to normal. The wait time may vary depending on your individual physiological characteristics.
10. Although such cases are rare, for those with an extremely weak pulse or irregular pulse, errors may result which prevent proper measurement. If abnormal variations are noticed, consult with your physician or trained healthcare professional.
11. This device is intended for adult use. While taking a measurement, you can stop the inflation or deflation process of the cuff at any time by pressing the POWER button.

12. If any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user is established.
13. Do not expose the device to extreme temperatures, humidity dust or direct sunlight as this may cause it to malfunction.
14. Please comply with the storage and operating conditions defined in 'Technical Specification' . Storing or using the device outside of the specified temperature and humidity range can affect measurement accuracy or the function of the device.
15. If the device was not stored within the minimum/maximum permissible storage conditions, a waiting period of at least 2 hours must be observed before using it under the specified operating conditions ( 'Technical Specification' ) or an ambient temperature of approx. 20 °C.

For Customer Service, the blood pressure monitor is calibrated when manufactured; it is recommended that the accuracy should be maintained and calibrated by manufacture triennially (every 3 years). To obtain the service please contact AViTA Corp. for the address of the repair location. Enclose with the Proof of Purchase. Include \$10.00 USD for the return shipping and handling. Accompany with a letter, with your name, address, phone number, and description of the specific problem or routine check-up. Pack the device carefully with bubble wraps (if there is) to prevent damages cause during transit. Due to possible losses in transit, it is recommended insuring the device with return receipt requested. If in any way of assistance of setting up, using, maintaining or to report unexpected operation/adverse events please contact manufacturer or local representative for further assistance.

 **CAUTION:**

- Beware of blood flow interference and resulting harmful injury to the patient caused by continuous cuff pressure.
- While using monitoring ME equipment and blood pressure monitor simultaneously being used on the same limb, pressurization of the cuff may cause monitoring ME equipment temporarily dysfunction.
- Reading can be affected by the measurement site, the position of the patient, exercise, or patient's physiologic condition.
- Automated sphygmomanometer can be affected by extremes of temperature, humidity and altitude.
- Frequent measurements can cause injury to the patient due to blood flow interference.
  
- Improper operation of automated sphygmomanometer may result in prolonged impairment of patient blood circulation, please read manual for proper use, if still not clear, contact professional health care personnel or local distributor for assistant.
- Retake the measurement if unexpected readings are obtained.
- Ensure that the cuff is not placed on an arm in which the arteries or veins are undergoing medical treatment, e.g., intravascular access or therapy, or an arteriovenous (AV) shunt.
  
- This product is suitable for use in the home healthcare environment.
- Keep this device out of the reach of children. Strangulation may result from baby or child entanglement in cables.
- Please keep this device away from pets, pests, and children.
- Preventing potential allergic reaction, please avoid the device in direct contact to patient's wound.

- Do not use the cuff on people who have undergone a mastectomy.
- Do not apply the cuff over a wound, as this can cause further injury.
- Do not apply the cuff other than the original manufacturer provided.
- Do not use in these cases (e.g. Device for use in an ambulance, helicopter or professional environment)
- Cuff pressure 0 – 300 mmHg
- Reduction rate:  $\leq 30\text{S}$  Refer to IEC 80601-2-30
- No modification of this equipment is allowed.
- High BMI health condition users may result varied.
- Users shall notify serious adverse event to the central competent authority or its commissioned agency, legal entity, or manufacture.

## Product Identification



- 1 Display
- 2 Cuff
- 3 Memory button (M key)
- 4 Start/Stop button (Power key)
- 5 Set button (S key)



## Description of LCD Display



Low battery indicator



WHO indicator



Date & Time



Irregular Heartbeat Symbol



AF Symbol



Memory Symbol



Memory Set



Heartbeat Symbol



Pulse rate



Inflation / Releasing air



Systolic Pressure




Diastolic Pressure



User Symbol

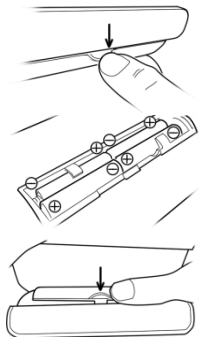
## Battery Installation

### Low battery warning:

It is necessary to replace the batteries when the Low Battery symbol "" appears on the display, or when the display does not turn on after the POWER key is pressed.

### Replacing the Battery:

1. Press down on latch and lift the cover on the bottom of the monitor.
2. Insert or replace 4x 1.5 V AAA batteries into the battery compartment, ensuring to match the indicated polarity symbols.  
Always use new batteries.
3. Replace the battery cover.



### CAUTION:

Risk of explosion if battery is replaced by an incorrect type. Dispose of used batteries according to the instructions.

### NOTE: Battery-operated

1. Please properly dispose of the batteries away from small children and heat. Avoid the children accidentally swallow the battery.
2. It is recommended to remove the batteries if the device will not be used for more than 1 month.
3. Batteries must be disposed of in accordance with local environmental and institutional policies.
4. It is recommended not to use rechargeable, unqualified or different spec battery may damage the device or cause circuit shortcut.

## Setting the Date and Time

It is necessary to set the date and time for the device every time batteries are initially installed or replaced.

To set the date and time, proceed as follows:

- While in power off mode, press and hold the “S” key for at least 3 seconds to enter Date and Time setting procedure and the Year value will begin to flash.
- Press the “S” key to advance the display to the desired year, press the “M” key to confirm the year.
- Next, the month value will blink. Repeat step 2 to set the month and date, then hours, then minutes.
- After setting the minutes, the device will automatically exit out of the date/time setting mode and shut off.

## Placement of the Pressure Sleeve

It is important to avoid smoking, eating, taking medication, alcohol consumption or physical activity 30 minutes prior to taking a reading. If for any reason you are unable to or should not use your left arm, please modify the instructions for cuff application to your right arm. Your physician can help you identify which arm is best for you to take measurements from.

1. Remove any constrictive clothing or jewelry that may interfere with cuff placement.
2. Be seated at a table or desk with your feet flat on the floor.

### **Note:**

Blood pressure naturally varies from one arm to the other; therefore, measure your blood pressure on the same arm to ensure comparability of the two readings.

3. Fully open the cuff by pulling the bottom of it. Then grip the clip of the device to widen the cuff, creating a clamp shape.
4. Insert your left arm into the created clamp shape, with the direction indicator align with your brachial artery.
5. The bottom edge of the cuff should be positioned approximately one inch above the elbow joint.
6. Reach underneath your left arm with your right hand, and pull the end of the cuff towards your body to tighten the cuff. Wrap and secure the cuff.
7. The cuff should fit comfortably, yet snugly around your arm. You should be able to insert one finger easily between your arm and the cuff.

**Note:**

If you are not comfortable with applying your cuff, please seek the assistance of another member of your household or work with your physician to practice the cuff application. Incorrectly applied cuffs may result in inaccurate readings.

## Measurement of Pulse Rate and Blood Pressure

1. Press the Power key to turn the power on. After full display is shown, the values for the last reading will appear on the display. If there is no measurement, the unit displays the value "0" .
2. AF mode selection will appear after the value of last reading. Press the Power key again to start the measurement.
3. After the self-test, the blood pressure monitor starts to measure. The cuff will automatically begin to inflate, with the display showing the increasing pressure in the cuff.
4. As the pressure increases, the indicator will increase upwards according to the pressure value on the display.
5. As the cuff inflates, the monitor automatically determines your ideal inflation level. This monitor detects your blood pressure and pulse rate during inflation.
6. To detect the heartbeat, the heartbeat symbol will appear and continuous flashes on the LCD display.
7. Your blood pressure measurement and pulse will display simultaneously on the screen.
8. The Hypertension Indicator will indicate your reading range on the display separately.
9. Press the Power key to turn the device off, conserving energy and battery life. The device will automatically shut-off after approximately 2 minutes of idling.

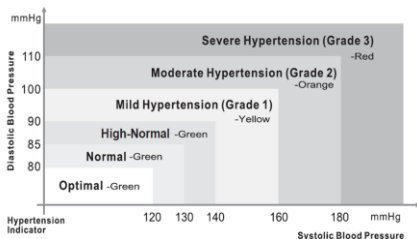
If there is an Error message after the measurement, please refer to Page 18 Error Codes Table.

### Note:


1. Users are recommended to use BPM67B daily at a fixed time, frequency of usage from 1 to 3 times can be suggested from physician or doctor' s advice.
2. Base on individual body condition the typical operation time takes approximately 30 seconds to 1 minute.

## Hypertension Indicator

This device features our unique Hypertension Indicator. The World Health Organization has established globally accepted standards for the assessment of high or low blood pressure readings. The below chart should be considered only as a guideline, always consult with your physician or health care professional to interpret your individual results.



## Irregular Heartbeat Detector

Your digital blood pressure monitor features an Irregular Heartbeat Detector. This feature allows users to accurately monitor blood pressure even if an irregular heartbeat should occur. When an irregular heartbeat is detected, the icon “” will appear on the display.

### Note:

Please consult with your physician or trained healthcare professional for further information regarding an irregular heartbeat if this symbol appears frequently.


### IMPORTANT INFORMATION:

This blood pressure monitor is not designed for use by people with arrhythmias nor for diagnosing or treating an arrhythmia problem. As a safeguard, we recommend that if you have arrhythmias such as atrial or ventricular premature beats and atrial fibrillation or any other special conditions you should check with your physician before using your blood pressure monitor.

## Atrial Fibrillation (AFIB)

This device is a blood pressure monitor that also analyses heart rate variability during measurement.

While in power off mode, press Power key into select AF mode, and use “S” key to select “ON” or “OFF” the AF mode. Confirm your selection with the Power key to turn the power on and start to measure your blood pressure.

This feature allows users to accurately monitor blood pressure even if atrial fibrillation should occur. When atrial fibrillation is detected, the “” icon will appear on the display.



Note:

1. Please consult with your physician or trained healthcare professional for further information regarding an atrial fibrillation if this symbol appears frequently.
2. Even if the AF symbol does not appear, there is still a possibility of Afib.

## Memory Function

### Recalling Measurements in Memory:

You can recall up to 100 measurements for each user, plus those average values stored measurements in memory to share with your physician or trained healthcare professional.

If you press the M key, the unit will first display the average of last 3 currently stored measurements. Continue to press the M key to successively view the average value of the morning measurements and evening measurements for the last 7 days will be displayed.

The duration of morning is AM5:00 – AM9:00,

The duration of evening is PM6:00 – PM8:00

If you press the M key again, the measurements will appear on the display from most current to oldest. And the memory number will appear on the display.

All results for a given measurement will display, including measurement results, pulse rate, Hypertension Indicator, Irregular Heartbeat alert, and date/time stamp.

When the number of readings exceeds 100, the oldest data will be replaced with the new record.

Press the Power key to turn the device off at any time when you check the stored measurements

### Clearing Measurements from Memory:

From the memory mode, press and hold down the M key until the display shows CLr. This indicates that all measurements have been erased.

## Data Transferring

This product is a arm type blood pressure monitor. Design without entering personal information. The device has a transmission function. During the transmission process, the measurement data is designed to be encrypted, and will not be tampered with unauthorized access nor retrieve user-related information. The firmware and software of the product have been pre-programmed in the production stage. Both the burned in microcontroller programmer and the software were designed encrypted; there is no concern about the software safety.

Bluetooth function requirement:

An Android device with Android version 4.3 or above and hardware support for Bluetooth 4.2.

An iOS device with iOS version 5 or above and hardware support for Bluetooth 4.2.

Note:

Please refer to the instruction manual of your smart phone for how to activate the Bluetooth function.

### Set Up Process

- (1) Please determine your mobile phone or computer has BLE4.2.
- (2) Turn on the device, when the Bluetooth icon shows on the device, it means device is under broadcast condition.
- (3) Please check the connecting condition from your mobile phone or computer.
- (4) Every measure reading will be transfer to your mobile phone or computer automatically.

You can send results in memory zone as following step:

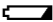
1. Press Memory recall button of blood pressure monitor to enter memory mode. The symbol AVG and M will appear on the display and the number of memory index is "A3" .
2. Turn on the Bluetooth function of device to search the blood pressure monitor. Then the blood pressure monitor will be searched and connected the Bluetooth device.
3. After the connection is completed, you can press memory recall button to display the memory result on LCD. Each result showed on LCD, the result will be transmitted by Bluetooth immediately.

## Two users with one device

This device is designed to be capable of recording the measurements of two people separately. To switch to the other user, please follow the steps below:

1. Press the "S" key while in power off mode and the user number will appear below the Date & Time on the display.
2. Change the user number by pressing the "S" key.
3. After selecting the user, press the power key to confirm and proceed to the measuring steps.
4. For checking the measurement records, press the "M" key after selecting the user.

## Error Codes

Err Code	Meaning	Corrective Action
Err 00	No pulse or detect pulses not enough.	Take off heavy clothes and retry again.
Err 01	The cuff is not fastened correctly, cuff pressure leakage or inflation too low or Overpressure protection.	The Arm cuff is not fastened properly. Re-apply the cuff, and take a measurement again.
Err 02	Inaccurate reading	Rest a while, relax and retry again.
Err 03	Inflation or deflation fail during the measurement	The Arm cuff is not fastened properly. Re-apply the cuff, and take a measurement again.
Err	Memory error.	Take off batteries to reboot the device, and then take another measurement. If the Err code persists after reinstalling the batteries, please contact the supplier or manufacturer.
	Low batteries	Replace all batteries with new ones.

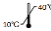
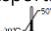
## Troubleshooting

Problem	Probable Cause	Recommended Action
Nothing appears in the display even when the power is turned on.	Batteries are drained.	Replace all batteries with new ones.
	Battery are not correctly aligned with terminals.	Reinsert batteries in the correct position.
Low Battery Symbol appears.	Batteries are drained.	Replace all batteries with new ones.
	In colder temperatures batteries have weaker electrical charges.	Warm up the batteries, or use the device in a warmer setting.
Device operation time is inconsistent.	Different battery brands have different life spans.	Use Alkaline batteries and replace all batteries at the same time with same brand batteries.
No reading after measurement.	Batteries are drained.	Replace all batteries with new ones.
Suspicious blood pressure results.	Perhaps the cuff was improperly positioned.	Adjust patient and Arm cuff to measure.
	Blood pressure naturally varies throughout the day.	Rest a while, relax and measure again.
Suspicious pulse rate results.	Bodily movement during device use.	Refrain from moving during measurement.
	Measurement shortly after exercise or exposure to the outdoors.	Do not take measurements after exercise or coming back from the outdoors.
Power switches off automatically.	System design.	Push the power button again, and then begin measure again.
During measuring, air re- inflates.	It could be a normal action if the user's blood pressure is higher than the initial pressure value, the device automatically pumps to a higher pressure by 40mmHg each time.	Relax, and try to take a measure again.
	The Arm cuff is not fastened properly.	Check that the Arm cuff is fastened properly and retake the measurement.

## Cleaning and Disinfecting

- Only use a soft, damp cloth to clean the monitor. Please do not use thinner, alcohol, detergents or solvents.
- The cuff can be cleaned carefully using a slightly damp cloth and mild soap solution. Do not completely immerse the cuff in water.
- It is recommended to clean and disinfect the cuff regularly or after each use, especially when used by several users, to prevent infection. The cuff should be disinfected, particularly on the inside, by wiping with a disinfectant. Use a disinfectant that is compatible with the cuff materials, e.g. 75 % ethanol or isopropyl alcohol.
- Keep the monitor in the appropriate carriage to protect it from external influences. And store it at appropriate condition.

## Technical Specification

- Measuring range :  
Blood Pressure : 40~255 mmHg  
Pulse Rate : 40~199 beats/min
- Calibration Accuracy:  
Blood Pressure :  $\pm 3$  mmHg  
Pulse rate :  $\pm 4\%$  of reading
- Operating environment :  
10°C~40°C   
15% to 85% relative humidity (non-condensing)  
700-1060 hPa ambient pressure
- Storage/ Transportation environment :  
-20 to 50 °C   
15% to 85% relative humidity (non-condensing)  
700-1060 hPa ambient pressure
- Power Source : 4 x 1.5 V AAA batteries
- Body Weight :  
approx. 250g (exclude batteries) +/- 5%
- Body Dimensions :  
approx. L 130 mm x W 94 mm x H 33.6 mm
- Cuff circumference (M/L Size) :  
approx. 22 ~ 42 cm (9" ~ 17" )
- Lifetime : 3 years
- Expected service life : 10,000 measurements

Note: After 3 years life time or 10,000 measurements, device material may experience degradation, measurement accuracy may vary.



















## EMC Tables








<b>BPM67B</b> is intended for use in the electromagnetic environment specified below. The customer or the user of <b>BPM67B</b> must make sure that it is used in such an environment.			
Guidance and manufacturer' s declaration - Electromagnetic emissions			
Phenomenon	Professional healthcare facility environment a)	HOME HEALTHCARE ENVIRONMENT a)	
Conducted and radiated RF MISSIONS	CISPR 11 Group 1 Class A (Not BLE Function) Group 2 Class A (With BLE Function)	CISPR 11 Group 1 Class B (Not BLE Function) Group 2 Class B (With BLE Function)	
Harmonic distortion	Not applicable b) (Note: Power by Battery or DC Input) Only the AC input needs to be tested		
Voltage fluctuations and flickering	Not applicable b) (Note: Power by Battery or DC Input) Only the AC input needs to be tested		
a) The equipment is suitable for use in Home Health Environments and Professional Health Care Environments limited to patient rooms and respiratory treatment facilities in hospital or clinics. The more restrictive acceptance limits of Group 1 Class B (CISPR 11) have been considered and applied. The equipment is suitable for use in the mentioned environments when directly connected to the Public Mains Network.			
b) The test is not applicable in this environment unless the ME EQUIPMENT and ME SYSTEM used will be connected to the PUBLIC MAINS NETWORK and the power input is otherwise within the scope of the Basic EMC standard.			
Guidance and manufacturer' s declaration - Electromagnetic immunity - Enclosure port			
Phenomenon	Basic EMC standard or test method	Immunity test levels	
		Professional healthcare facility environment	HOME HEALTHCARE ENVIRONMENT
ELECTROSTATIC DISCHARGE	IEC 61000-4-2	± 8kV contact ± 2 kV, ±4kV ±, ±8 kV, ±15 kV air	
Radiated RF EM fields	IEC 61000-4-3	a)	10 V/m b) 80MHz - 2.7 GHz 80% AM at 1kHz

Proximity fields from RF wireless communications equipment	IEC 61000-4-3	<p>COMPLIANT</p> <p>NOTE: Further information about distances to be maintained between portable and mobile RF communications equipment (transmitters) and the <b>BPM67B</b> can be requested from supplier using the contact information provided in this manual. However, it is advisable to keep the equipment at an adequate distance of, at least, 0.5 m from mobile phones or other RF communications transmitters to minimise possible interference.</p>
RATED power frequency magnetic fields.	IEC 61000-4-8	<p>30 A/m c) 50 Hz or 60 Hz</p>
<p>a) The equipment is suitable for use in Home Health Environments and Professional Health Care Environments limited to patient rooms and respiratory treatment facilities in hospital or clinics. The more restrictive IMMUNITY acceptance limits have been considered and applied.</p> <p>b) Before modulation is applied.</p> <p>c) This test level assumes a minimum distance of at least 15 cm between the ME EQUIPMENT or ME SYSTEM and sources of power frequency magnetic fields.</p>		

## Explanation of Symbols

Symbol	Definition
	The CE marking with the Registration Number of the Notified Body. This denotes the compliance of Regulation (EU) 2017/745
	Medical Device
	Manufacturer
	Authorized representative in the European Community
	Date of manufacture (YYYY-MM-DD or YYYY-MM)
	Batch code (YYMMWWWW)
	Serial number (YYMWWWXXXXX)
	Keep dry
	Temperature limit

	Humidity limitation
	Atmospheric pressure limitation
	Caution
	Consult the instruction for use
	Disposal information: Should you wish to dispose of the article, do so in accordance with current regulations. Details are available from your local authority. WEEE 2012/19/EU Directives
RoHS	This product fulfilling the requirements of the RoHS Directive 2011/65/EU.
REACH	This product fulfilling the requirements of the REACH Directive EC 1907/2006 and its amendments, do not contain Substances of Very High Concern in concentration above the limit of 0.1 %. No substance(s) is/are present in the parts of the product above the concentration of 0.1 % weight by weight.
	Stand-by
	Device classification type BF

<p><b>IP 22</b></p>	<p>This product meets the basic safety and essential performance requirements indicated in the IP22 conditioning test (protection against solid foreign objects of 12.5mm <math>\emptyset</math> and greater and against vertically falling water drops when enclosure tilted up to 15°)</p>
	<p>The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.</p>
	<p>Importer</p>
	<p>Distributor</p>
	<p>Model Number</p>
	<p>Country of Manufacturer</p>
	<p>Unique Device Identifier</p>
	<p>Keep away from sunlight</p>

Electronic IFU available at <http://www.avita.com.tw>

BPM67BP-22584AV  
2022-12-22