

TOUCH FREE Infrared Thermometer

NT17B



EN



CE
2797



IP 22 RoHS REACH



MDSS GmbH
Schiffgraben 41
30175 Hannover Germany



AVITA Corporation
9F., No. 78, Sec. 1, Kwang Fu Rd.,
San Chung Dist., New Taipei City 24158, Taiwan, China



China, AVITA(WUJIANG)



Importer



Distributor

INSTRUCTION MANUAL

Please read this instruction manual carefully
before using your forehead thermometer

Contents

Introduction	1
Important Information before Use	2
Product Identification	4
Description of LCD Display	5
Battery Installation	6
Switching Between Fahrenheit or Celsius	7
Switching between 4 Kinds of Measuring Mode.....	7
Tips for Measuring Human Temperature	8
Measuring Human Temperature.....	9
Measuring Object/Liquid Temperature	10
Memory Function	11
Cleaning and Disinfecting.....	12
Error Codes	13
Troubleshooting.....	14
Operating the Bluetooth function.....	16
Technical Specification	17
EMC Tables.....	18
Explanation of Symbols.....	20

Introduction

Intended Use:

Touch Free Infrared Thermometer is intended for home use and the intermittent measurement of human body forehead temperature in people of all ages.

Features:

1. Multi functions
Human Body/ Object/ Night Mode Option/ Fever Alarm/ Backlight
2. Night Mode Option
Switch to Night mode can mute the sound to reduce the interference of buzzer during your baby sleeping.
3. LED Light for Fever Alarm
This device has the LED light to remind user their readings are over 38°C and plus the beep sound only in Human Body Mode.
4. 30 sets Memory
5. Illuminated Backlight Display
6. °C/°F Switch-able Function
7. One-second Reading
8. Auto power off for power saving
9. Low-battery Indicator
Indicate for battery condition and measuring range.
10. Large LCD Display



Contraindications:

1. Scarred or compromised tissues on measurement site.
2. Patient' s skin' s on certain medication of measurement site.

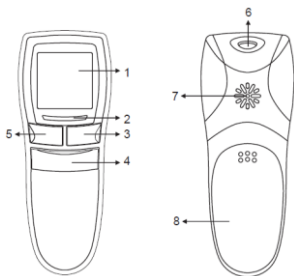
Important Information before Use

When using this product, please be sure to follow all the notes listed below. Any action against these notices may cause injury or affect the accuracy.

1. Do not disassemble, repair, or remodel the thermometer.
2. Be sure to clean the thermometer lens each time after usage.
3. Avoid direct finger contact with the lens.
4. No modification of this equipment is allowed.
5. It is recommended that user may take 3 temperatures. If they are different, use the highest reading.
6. Do not expose the thermometer to extreme temperature, very high humidity, or direct sunlight.
7. Avoid extreme shock or dropping the device.
8. Before the measurement, patients and thermometer should stay in steady state room condition for at least 30 minutes.
9. Avoid measuring temperature in 30 minutes after exercise, bathing, or returning from outdoor.
10. To protect the environment, dispose of empty batteries at appropriate collection sites according to national or local regulations.
11. Please use the thermometer solely for its intended purpose.
12. There are no absolute body temperature standards. Keep reliable records of your personal temperature to serve as a reference for judging a fever.
13. Under any circumstances, the temperature taking result is ONLY for reference. Before taking any medical action, please consult your physician.
14. The thermometer is calibrated at the time of manufacture. If at any time you question the accuracy of temperature measurements, please contact an authorized distributor for sending calibration services, additional cost may apply for return shipping fee.

15.  CAUTION: Please keep this device away from pets, pests, and children.
16.  CAUTION: Please keep away from children and pets, because small parts may be inhaled or swallowed.
17. If any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user is established.
18. If you would require any assistance during setting up, using or maintaining, please contact the manufacturer.
19. Avoid using the device near an operating cell phone or microwave oven.

Product Identification



- 1 LCD Display
- 2 LED Light
- 3 Mode Button
- 4 Measure Button
- 5 Memory Button
- 6 Measurement Sensor
- 7 Sound Hole
- 8 Battery Compartment Lid

Description of LCD Display



Measuring sequence in process



Night Mode (Mute Mode)



Human Mode



Object Mode



Degree of measurement



Celsius Scale



Fahrenheit Scale



Low battery indicator



Frowning Face
(temperature is equal or higher than 38°C)





Smiling Face
(temperature is lower than 38°C)



Memory Symbol and sets of records

Battery Installation

Low battery warning:

When the battery power becomes low, the low battery symbol “” will appear on the display. The thermometer can still be used during this time, but the batteries should be replaced as soon as possible. If the batteries run out completely, “Lo” will be displayed along with the low battery symbol “”. In this case, the batteries will need to be replaced before using the thermometer again.

CAUTION:

1. Risk of explosion if battery is replaced by an incorrect type.
2. It is recommended not to use rechargeable, unqualified or different spec of batteries as it may damage the device or cause circuit shortcut.

Replacing the Battery:

1. Gently slide the battery cover back.
2. Carefully remove the old batteries and properly discard.
3. Insert new batteries (Two 1.5V alkaline AAA(LR03) Size) according to the proper polarity.
4. Slide the battery cover back on.

NOTE: Battery-operated

1. Please properly dispose of the batteries away from small children and heat.
2. It is recommended to remove the batteries if the unit will not be used for an extended period of time.
3. Batteries must be disposed of in accordance with local environmental and institutional policies.
4. Dispose of used batteries in accordance with the applicable legal regulations. Never dispose of batteries in the normal household waste.

Switching Between Fahrenheit or Celsius

Your thermometer can display results in either degrees Celsius (°C) or degrees Fahrenheit (°F). With combination of Memory and Mode button are able to change the mode to either °C or °F with a beep sound.

Switching between 4 Kinds of Measuring Mode

1. Under power on status, you can press the **Mode Button** to switch different measuring mode. There are 4 kinds of mode which including Human, Object, Human/Night, and Object/Night Mode (in order).
2. The beep sounds will be off (mute) when your choice in Human/Night Mode, and Object/Night Mode, and the Moon symbol "☾" will appears on the LCD in both Night Modes.



NOTE:

Each press will come with a beep sound to ensure the setting is activated. (Except both Night Modes)

Tips for Measuring Human Temperature

Bear in mind that the thermometer needs to have been in the room in which the measurement is taken for at least 30 minutes before use.

NOTE:

- Attempting to take temperature readings from sites on the body other than the forehead may produce inaccurate results.
- The patient should remain still while the reading is being taken.
- Readings taken while asleep should not be compared directly to readings taken while awake, as body temperature while asleep is typically lower.
- Do not take body temperature readings within 30 minutes of being outdoors, exercising or bathing.

Measuring Human Temperature

Taking a Measurement

1. Press the Power Button to power on the thermometer. The unit will run a self-test and the LCD will briefly display all of its symbols during this time. When the device is ready, "--" appear on the screen, and the thermometer is ready to take a measurement.
2. Select the desired mode by pressing and releasing the Human/Object Button.


Note:

- When taking a patient's temperature, ensure that the thermometer is in Human Mode; the Human symbol "👤" will appear on the display.
 - If the patient's skin is covered with hair, sweat, or dirt, clean the area and wait 10 minutes before taking a measurement.
 - Ensure that the thermometer is held firmly during measurement and that the patient does not move until the measurement is complete. Movement can impact the measurement.
3. Position the thermometer under 3cm (around 1 inch) from the center of the patient's forehead with the sensor aimed between the eyebrows.
 4. Press and release the Measure Button. As you hear a short beep means this temperature reading has been completed and accompany with a back-light.
 5. If the temperature measurement is below 38°C, a "Smiling Face" will appear next to the reading. If the reading is 38°C or above, a "Frowning Face" will display and the LED light up to alert.
 6. After about 30 seconds after use, the thermometer will automatically beep and shut off.

Fever Alarm:

Only in Human Mode and Human/Night Mode.




Measuring Object/Liquid Temperature

1. Press the "Measure Button" to turn the thermometer on. All symbols on the display will momentarily appear.
2. Ensure that the thermometer is in Object Mode; the Object symbol  will be on the display. To alternate between modes press and release the Human/Object Mode Button until you see the desired measurement symbol on the display.
3. Position the thermometer under 3cm (around 1 inches) from the object.
4. Release the "Measure Button" and the temperature reading will be displayed.
5. After use, the thermometer will automatically beep and shut off after idling.

Memory Function

Memory Recall:

You can recall up to 30 measurements currently stored in memory to share with your physician or trained healthcare professional.

1. When the device is on, press once briefly on the Memory Button, then pass it again to show the last measurement accompanied by “” symbol.
2. The “” symbol or “” symbol will appear with each measurement stored in memory to indicate whether a person or object temperature was taken.
3. Each press of the same button recalls a previous measurement.

Memory Deletion:

1. Under power on status, you may hold press the “Memory button” to delete all the readings. The LCD shows “- -” with one beep sound to indicate that all memories are cleared.
2. Remove batteries, memory will be cleared.

NOTE:

All the readings will be cleared no matter record in Human Mode or Object Mode.

Cleaning and Disinfecting

For home use device disinfection, 75% alcohol (available in the pharmacy) can be used.

- **Measurement Sensor**

Clean the measurement sensor with an alcohol swab before and after each measurement.

- **Thermometer:**

Use a soft, dry cloth to clean thermometer body. Never use abrasive cleaning agents, thinners or benzene for cleaning. Do not scratch the surface of the probe lens or the display. Do not expose the thermometer to extreme temperatures, humidity, direct sunlight, or shock.

Error Codes

When a malfunction or incorrect temperature measurement occurs, an error message will appear as described below.

Display	Cause	Solution
Hi	The temperature measured is higher than measurement range. Human mode: 43°C (109.4°F) Object mode: 100°C (212.0°F)	Operate the thermometer only between the specified temperature ranges. If necessary, clean the sensor tip. In the event of a repeated error message, contact your retailer or Customer Services.
Lo	The temperature measured is lower than measurement range. Human mode: 34°C (93.2°F) Object mode: 0°C (32.0°F)	
Err	The operating temperature is not in the range: 15°C~40°C (59°F~104°F)	Operate the thermometer only between the specified temperature ranges.

Troubleshooting

Trouble	Probable Case	Recommended Action
Failed to power on.	Batteries are drained.	Replace with a new battery
	Batteries are not correctly aligned with terminals.	Reinsert batteries correctly.
	The thermometer is damaged.	Contact your retailer or Customer Services.
Low battery symbol appears.	Low battery.	Replace with a new battery as soon as practical
	In colder temperature batteries have weaker electrical charges.	Warm up the batteries or use the device in a warmer setting.
The measurement is abnormal or if there is any doubt on the measured result	Incorrect measure position.	Reposition the probe point center of the forehead.
	The measuring distance is too far.	Maintain a distance of 3 cm.
	The measuring is hampered by hair, perspiration, oil, sweat, makeup, etc.	Remove obstructions from forehead.
	The lens of probe is dirty.	Clean the lens according to <i>cleaning and disinfecting</i> .
	You have just come	Stay in steady state

	from extreme temperatures environment.	room at least for 30 minutes.
	The thermometer has been stored in a cold room or exposed to direct sunlight.	Put the thermometer in steady state room at least for 30 minutes.
	Measurement was taken in surface temperature mode.	Switch mode to human mode.
Power switches off automatically.	System design.	Push the power button again.
	Batteries are drained.	Replace with a new battery
	The thermometer is damaged.	Contact your retailer or Customer Services.

Operating the Bluetooth function

This product is an **TOUCH FREE Infrared Thermometer with Bluetooth function. design** without entering personal information. If the device has a transmission function, the transmission measurement data is designed to be encrypted and transmitted, and will not be tampered with or retrieve user-related information during the transmission process. The firmware and software of the product have been programmed in the production process, and the programming interface is different from the data transmission interface. When programming to the microcontroller, use an encrypted programmer, so there is no need to worry about the software being tampered with during transmission.

Bluetooth function requirement:

- An Android device with Android version 4.3 or above and hardware support for Bluetooth 4.2.
- An iOS device with iOS version 5 or above and hardware support for Bluetooth 4.2.

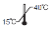

How to activate the Bluetooth function:

Please refer to the instruction manual of your mobile phone or computer for how to activate the Bluetooth function.

Set Up Process

- (1) Please determine your mobile phone or computer has BLE4.2.
- (2) Turn on the **NT17B device**, when the Bluetooth icon shows on the device , it means **NT17B device** is under broadcast condition.
- (3) Please check the connecting condition from your mobile phone or computer. The device name should be “**NT17B**” .
- (4) Every measure reading will be transfer to your mobile phone or computer automatically.

Technical Specification

- **Measuring range :**
Human Body : 34°C~43°C (93.2°F~109.4°F)
Object : 0°C~100°C (32.0°F~212.0°F)
- **Calibration Accuracy:**
Human Body : $\pm 0.2^{\circ}\text{C}(\pm 0.4^{\circ}\text{F})$: from 34°C~43°C (93.2°F~109.4°F)
Object : $< 40^{\circ}\text{C} \pm 2^{\circ}\text{C}; \geq 40^{\circ}\text{C} \pm 5\%$
- **Display resolution : 0.1°C**
- **Operating environment :**
15°C~40°C (59°F ~ 104°F) 
with relative humidity 15% to 95% (non condensing)
Atmospheric pressure: 700~1060 hPa
- **Storage/ Transportation environment :**
-25 to 55 °C(-13 to 131 °F) 
with relative humidity 15% to 95% (non condensing)
- **Atmospheric pressure: 700~1060 hPa**
- **Power supply : 2 x 1.5V AAA size alkaline batteries**
- **Weight : approx. 80g (with batteries)**
- **Memory: 30 sets**
- **Dimensions :**
approx. 128.5mm×48.83mm×38.85mm (L×W×H)
- **Operation Distance : 3 cm**
- **Frequency : 2402~2480GHz**
- **Output power range : $\leq 4\text{dBm}$**
- **Reference body site: Ear temperature**
- **Shelf Life : 3 years**
- **Battery Operation Time:**
about 90 minutes (Depends on different brand of batteries.)
- **For Customer Service**










To obtain further service please contact AViTA Corp. for the address of the repair location. Enclose the Proof of Purchase. Include \$10.00 USD for the return shipping and handling. Include a letter, with your name, address, phone number, and description of the specific problem. Pack the product carefully to prevent damage in transit. Because of possible loss in transit, we recommend insuring the product with return receipt requested.








EMC Tables




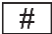



<p>NT17B is intended for use in the electromagnetic environment specified below. The customer or the user of NT17B must make sure that it is used in such an environment.</p>			
<p>Guidance and manufacturer's declaration - Electromagnetic emissions</p>			
Phenomenon	Professional healthcare facility environment a)	HOME HEALTHCARE ENVIRONMENT a)	
Conducted and radiated RF MISSIONS	CISPR 11 Group 1 Class A (Not BLE Function) Group 2 Class A (With BLE Function)	CISPR 11 Group 1 Class B (Not BLE Function) Group 2 Class B (With BLE Function)	
Harmonic distortion	Not applicable (Note: Power by Battery or DC Input) Only the AC input needs to be tested		
Voltage fluctuations and flickering	Not applicable (Note: Power by Battery or DC Input) Only the AC input needs to be tested		
<p>a) The equipment is suitable for use in Home Health Environments and Professional Health Care Environments limited to patient rooms and respiratory treatment facilities in hospital or clinics. The more restrictive acceptance limits of Group 1 Class B (CISPR 11) have been considered and applied. The equipment is suitable for use in the mentioned environments when directly connected to the Public Mains Network.</p> <p>b) The test is not applicable in this environment unless the ME EQUIPMENT and ME SYSTEM used will be connected to the PUBLIC MAINS NETWORK and the power input is otherwise within the scope of the Basic EMC standard.</p>			
<p>Guidance and manufacturer's declaration - Electromagnetic immunity - Enclosure port</p>			
Phenomenon	Basic EMC standard or test method	Immunity test levels	
		Professional healthcare facility environment	HOME HEALTHCARE ENVIRONMENT
ELECTROSTATIC DISCHARGE	IEC 61000-4-2	± 8kV contact ± 2 kV, ±4kV ±, ±8 kV, ±15 kV air	
Radiated RF EM fields	IEC 61000-4-3	a)	10 V/m b) 80MHz - 2.7 GHz 80% AM at 1kHz

Proximity fields from RF wireless communications equipment	IEC 61000-4-3	<p>COMPLIANT</p> <p>NOTE: Further information about distances to be maintained between portable and mobile RF communications equipment (transmitters) and the NT17B can be requested from supplier using the contact information provided in this manual. However, it is advisable to keep the equipment at an adequate distance of, at least, 0.5 m from mobile phones or other RF communications transmitters to minimise possible interference.</p>
RATED power frequency magnetic fields.	IEC 61000-4-8	<p>30 A/m c) 50 Hz or 60 Hz</p>
<p>a) The equipment is suitable for use in Home Health Environments and Professional Health Care Environments limited to patient rooms and respiratory treatment facilities in hospital or clinics. The more restrictive IMMUNITY acceptance limits have been considered and applied.</p> <p>b) Before modulation is applied.</p> <p>c) This test level assumes a minimum distance of at least 15 cm between the ME EQUIPMENT or ME SYSTEM and sources of power frequency magnetic fields.</p>		

Explanation of Symbols

	The CE marking with the Registration Number of the Notified Body. This denotes the compliance of Regulation (EU) 2017/745
	Medical Device
	Manufacturer
	Authorized representative in the European Community
	Date of manufacture (YYMMDD or YYYY-MM)
	Batch code (YYMMWWWW)
	Serial number (YYMWWWWXXXXX)
	Keep dry
	Temperature limit

	Humidity limitation
	Atmospheric pressure limitation
	Caution
	Consult the instruction for use
	Disposal information: Should you wish to dispose of the article, do so in accordance with current regulations. Details are available from your local authority. WEEE 2012/19/EU Directives
RoHS	This product fulfilling the requirements of the RoHS Directive 2011/65/EU.
REACH	This product fulfilling the requirements of the REACH Directive EC 1907/2006 and its amendments, do not contain Substances of Very High Concern in concentration above the limit of 0.1 %. No substance(s) is/are present in the parts of the product above the concentration of 0.1 % weight by weight.
	Stand-by
	Device classification type BF

<p>IP 22</p>	<p>This product meets the basic safety and essential performance requirements indicated in the IP22 conditioning test (protection against solid foreign objects of 12.5mm \varnothing and greater and against vertically falling water drops when enclosure tilted up to 15°)</p>																
	<p>The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.</p>																
	<p>Importer</p>																
	<p>Distributor</p>																
	<p>Model Number</p>																
	<p>Country of Manufacturer</p>																
	<p>Unique Device Identifier</p> <table border="1" data-bbox="262 826 967 1065"> <thead> <tr> <th data-bbox="262 826 335 882">AI</th> <th data-bbox="335 826 672 882">Data Element</th> <th data-bbox="672 826 967 882">Data Type/Format</th> </tr> </thead> <tbody> <tr> <td data-bbox="262 882 335 926">(01)</td> <td data-bbox="335 882 672 926">Device Identifier(GTIN)</td> <td data-bbox="672 882 967 926">Numeric</td> </tr> <tr> <td data-bbox="262 926 335 971">(11)</td> <td data-bbox="335 926 672 971">Manufacturing Date</td> <td data-bbox="672 926 967 971">Numeric: YYMMDD</td> </tr> <tr> <td data-bbox="262 971 335 1015">(10)</td> <td data-bbox="335 971 672 1015">Lot Number</td> <td data-bbox="672 971 967 1015">Alphanumeric</td> </tr> <tr> <td data-bbox="262 1015 335 1065">(21)</td> <td data-bbox="335 1015 672 1065">Serial Number</td> <td data-bbox="672 1015 967 1065">Alphanumeric</td> </tr> </tbody> </table>		AI	Data Element	Data Type/Format	(01)	Device Identifier(GTIN)	Numeric	(11)	Manufacturing Date	Numeric: YYMMDD	(10)	Lot Number	Alphanumeric	(21)	Serial Number	Alphanumeric
AI	Data Element	Data Type/Format															
(01)	Device Identifier(GTIN)	Numeric															
(11)	Manufacturing Date	Numeric: YYMMDD															
(10)	Lot Number	Alphanumeric															
(21)	Serial Number	Alphanumeric															
	<p>Keep away from sunlight</p>																

Electronic IFU available at <http://www.avita.com.tw>

NT17BP-22428AV

2023-05-30